

ADAPT

Transtibial / Transfemoral Liners

The ADAPT liner is a silicone gel liner with embedded mineral oil. The soft gel with flexible fabric allows snug fit of various contours without restriction. The liner conforms very nicely while providing comfort. It is suitable for moderate to high activity patients. The matrix which is built into the liner prevents distal elongation of the liner for locking application. This product has a six-month warranty.

Transtibial (LOCKING)



P/N	Description - LOCKING
3012L-18	ADAPT Locking Liner, 18
3012L-20	ADAPT Locking Liner, 20
3012L-22	ADAPT Locking Liner, 22
3012L-23.5	ADAPT Locking Liner, 23.5
3012L-25	ADAPT Locking Liner, 25
3012L-26.5	ADAPT Locking Liner, 26.5
3012L-28	ADAPT Locking Liner, 28
3012L-30	ADAPT Locking Liner, 30
3012L-32	ADAPT Locking Liner, 32
3012L-34	ADAPT Locking Liner, 34
3012L-36	ADAPT Locking Liner, 36



Transtibial (CUSHION)


P/N	Description - CUSHION
3012C-18	ADAPT Cushion Liner, 18
3012C-20	ADAPT Cushion Liner, 20
3012C-22	ADAPT Cushion Liner, 22
3012C-23.5	ADAPT Cushion Liner, 23.5
3012C-25	ADAPT Cushion Liner, 25
3012C-26.5	ADAPT Cushion Liner, 26.5
3012C-28	ADAPT Cushion Liner, 28
3012C-30	ADAPT Cushion Liner, 30
3012C-32	ADAPT Cushion Liner, 32
3012C-34	ADAPT Cushion Liner, 34
3012C-36	ADAPT Cushion Liner, 36

Transfemoral (LOCKING ONLY)




P/N	Description - LOCKING
3092L-30	ADAPT Locking TF Liner, 30
3092L-32	ADAPT Locking TF Liner, 32
3092L-34	ADAPT Locking TF Liner, 34
3092L-36	ADAPT Locking TF Liner, 36
3092L-38	ADAPT Locking TF Liner, 38
3092L-40	ADAPT Locking TF Liner, 40
3092L-45	ADAPT Locking TF Liner, 45


HOW TO MEASURE Transtibial Liner

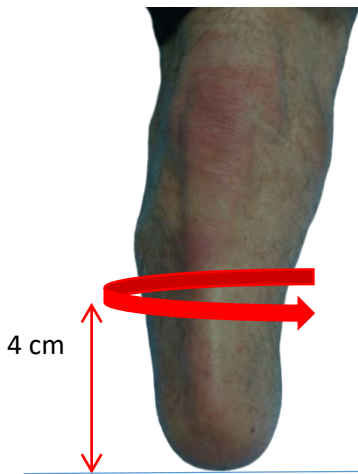


1. Measure the circumference of the residual limb at 4cm from the distal end
2. Choose the liner size, which corresponds to the measurement or one smaller.
3. Never choose a liner size, which is larger than the measurement

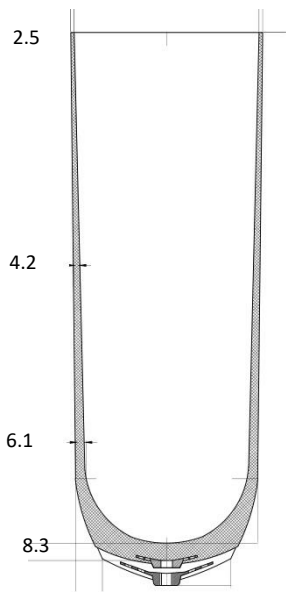


ADAPT LOCKNG

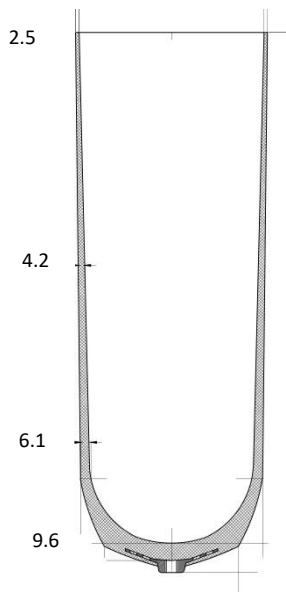

PARASIL LOCKING

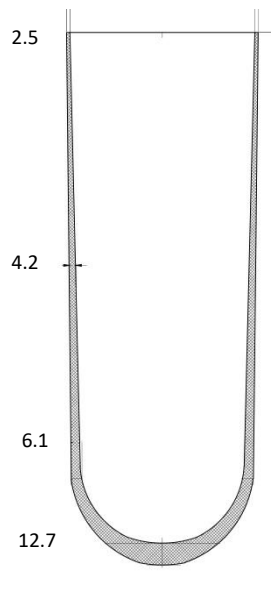

ADAPT/PARASIL
CUSHION

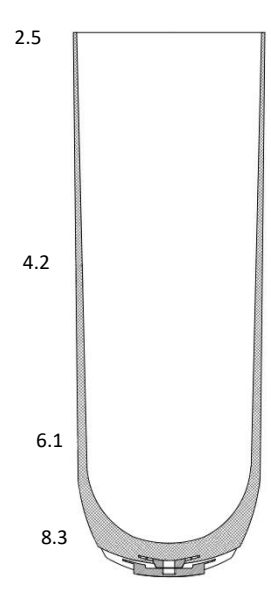


PREMIER/HYDRO
MAGLINER

BK LINERS


WALL AVERAGES

www.stngco.com


HOW TO MEASURE Transfemoral Liner

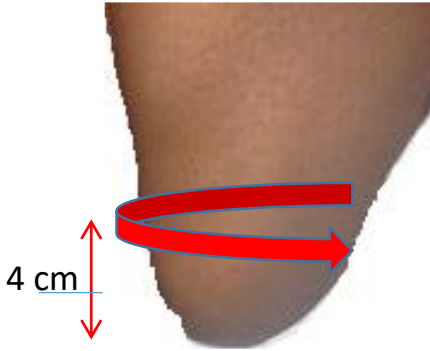


1. Measure the circumference of the residual limb at 4cm from the distal end
2. Choose the liner size, which corresponds to the measurement or one smaller.
3. Never choose a liner size, which is larger than the measurement

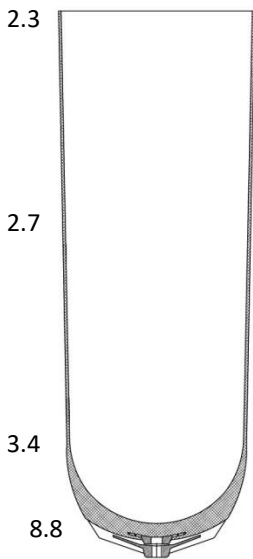


ADAPT LOCKING

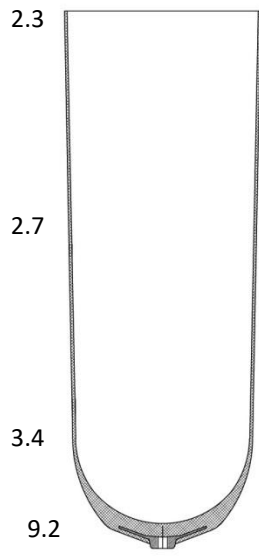

PARASIL LOCKING

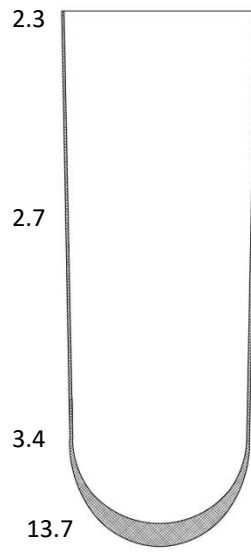

PARASIL
CUSHION

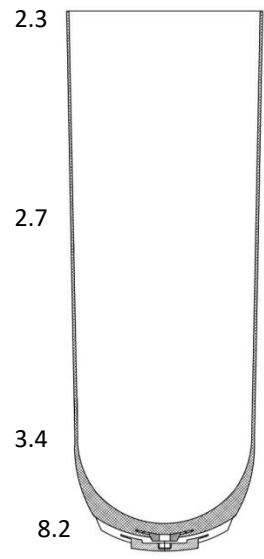

PREMIER/HYDRO
MAGLINER

AK LINERS

WALL AVERAGES

www.stngco.com